


# Obstruction Removal Kit for ST Nebulizers

On occasion, the ST nebulizers may become clogged. This kit is specifically designed for PFA-ST and PolyPro-ST nebulizers.

- Includes a guide tool and 10 pieces of fused silica.
- 0-3 ml/min sample flow rate with peristaltic pumping.
- Easy to use:
  - tighten guide onto back of ST nebulizer
  - slide fused silica through
- Silica refills also available.

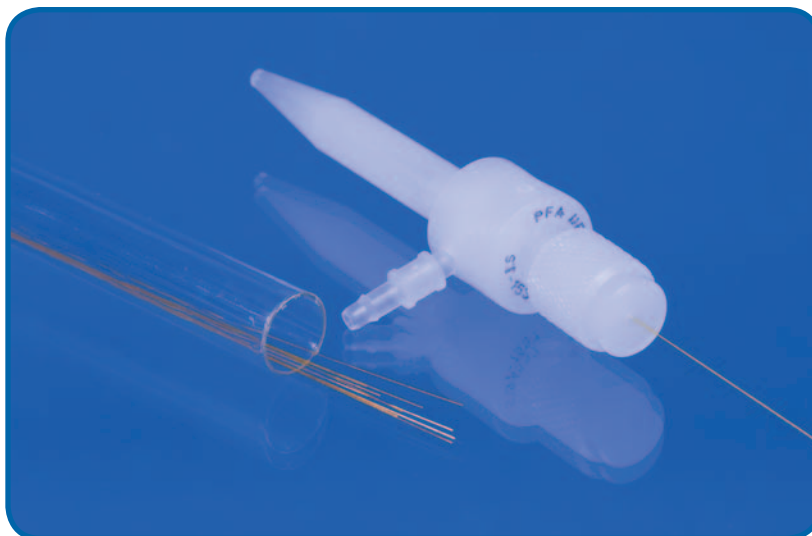

Obstruction Removal Kit shown with PFA-ST Nebulizer (not included)

Ordering Information	
<b>ES-2046-0001</b>	Fused Silica Obstruction Removal Kit for ST-type Nebulizers. Includes one guide, 10 pieces of fused silica and instructions.
<b>ES-2046-0002</b>	Fused Silica Refill for Obstruction Removal Kit For use with ST-type Nebulizers only. Package of 10.

